



Instructions for use

PAX Roll-Up Stretcher London Ambulance Service





PAX Roll-Up Stretcher London Ambulance Service



Product: PAX Roll-Up Stretcher
London Ambulance Service

Item Number: 268935803



Manufacturer: X-CEN-TEK GmbH & Co. KG
Westerburger Weg 30
D-26203 Wardenburg / Germany
Phone: +49 (0) 44 07 - 7 14 76 0
Fax: +49 (0) 44 07 - 7 14 76 99
Email: info@pax-bags.de
Website: www.pax-bags.com

Status: 20.07.2021
(Date of last revision)

Appropriate use / purpose

The PAX Roll-Up Stretcher London Ambulance Service is used to transport a patient lying down.

With the stretcher the patient can be transported by persons carrying it.

The patient can also be transported in a rescue vehicle. Max carry load: 200 kg

**To use the PAX Roll-Up Stretcher London Ambulance Service,
the rescue personnel must be specially trained in the use of the product.**

Attention

**Before use, training in the use
of this product
is essential.**

Liability

Before using the PAX Roll-Up Stretcher London Ambulance Service for the first time, it is essential that you read and understand the instructions carefully.

The PAX Roll-Up Stretcher London Ambulance Service should only be used by people in good physical and mental condition.

They must be trained in the safe use of the PAX Roll-Up Stretcher London Ambulance Service and have the necessary expertise to use the product or be under the supervision of such a person.

Attention

**Activities in which this equipment
is used, are inherently
dangerous.
All actions and decisions are the
responsibility of the user.**

Before using this equipment:

- The instructions for use must be read and understood completely,
- the user must be trained in the correct use of this equipment.
- The user must familiarise themselves with the equipment and understand the possibilities and restrictions.
- The user understands the risks associated with its use and accepts these risks.

**Failure to heed any of these warnings
can lead to serious injuries or even death!**




Instructions for use PAX Roll-Up Stretcher London Ambulance Service

Table of contents

Chapter	Page	Chapter	Page
General information	2	5) Information on the label	11
1) Component names	4	6) Compatibility	11
2) Before use / after use	6	7) Safety requirements	11
3) Packaging	8	8) Inspection and service life	12
A) Preparing the vacuum support	9	9) Maintenance and storage	12
B) Patient positioning	9	10) Cleaning	12
C) Preparation for transport	10	11) Modification, repair	12
4) General information	11	12) Control card	13



1) Component Names

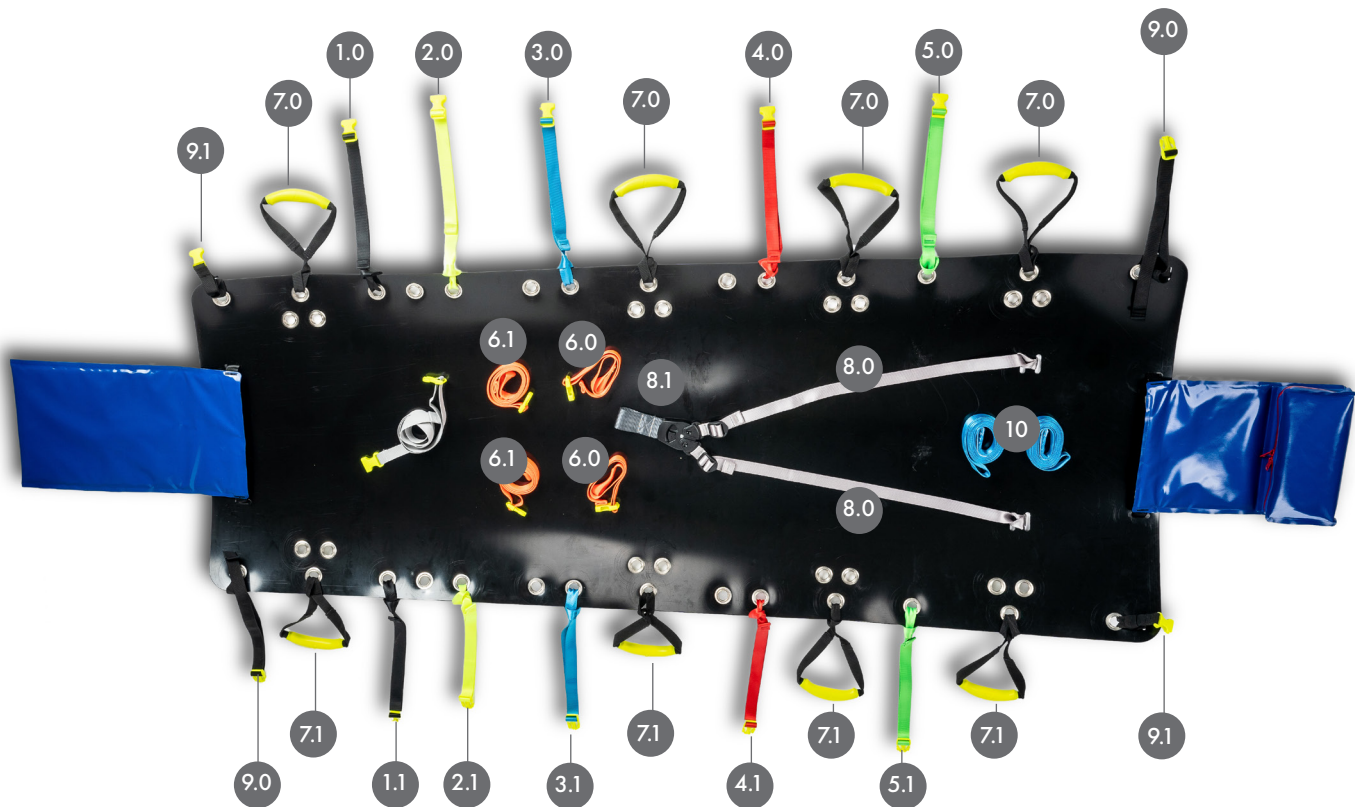

- | | | |
|------------------------|------------------------|----------------------------------|
| 1.0 Belt straps black | 4.0 Belt straps red | 7.0 Handles 1 |
| 1.1 Belt straps black | 4.1 Belt straps red | 7.1 Handles 2 |
| 2.0 Belt straps yellow | 5.0 Belt straps green | 8.0 Belt straps grey |
| 2.1 Belt straps yellow | 5.1 Belt straps green | 8.1 Belt straps grey |
| 3.0 Belt straps blue | 6.0 Belt straps orange | 9.0 Belt straps black |
| 3.1 Belt straps blue | 6.1 Belt straps orange | 9.1 Belt straps black |
| | | 10 PAX Sling light blue (150 cm) |



1) Component Names

item number: 270324200

**1 set PAX Roll-Up-Stretcher LAS
- fixation belt set**

The set includes 5 different coloured
25mm belt straps with high-viz yellow buckle
(adjustable on both sides)

Belt straps colours: blue, yellow, red, black, green

Belt straps lengths:

1.0	black	80 cm
1.1	black	80 cm
2.0	yellow	80 cm
2.1	yellow	80 cm
3.0	blue	80 cm
3.1	blue	80 cm
4.0	red	80 cm
4.1	red	80 cm
5.0	green	80 cm
5.1	green	80 cm

item number: 270334209

**1 set PAX Roll-Up-Stretcher LAS
- leg fixation belt set**

1 set consists of 2 complete belts
25mm belt straps with high-viz yellow buckle
(adjustable on both sides)

Belt straps colours: orange

Belt straps lengths:

6.0	orange	70 cm
6.1	orange	180 cm

item number: 162420008

**4 sets PAX Vacuum mattress handle set
(2 handles)**

1 set includes 2 handles

7.0	Handle 1
7.1	Handle 2

item number: 270374310

**1 piece PAX Spare part PAX Roll-Up-Stretcher LAS
- chest Fixtion**

25mm belt strap and 40 mm belt strap grey
with click buckle

Belt straps lengths:

8.0	grey	25 mm	89 cm	2 pieces
8.1	grey	40 mm	10 cm	1 piece

item number: 270381103

**1 set PAX Spare part PAX Roll-Up-Stretcher LAS
- mounting belt set**

25mm belt strap black with
standard high-viz yellow buckle.
1 set includes 2 belt straps

Belt straps lengths:

9.0	black	80 cm
9.1	black	20 cm

item number: 141034202


1 piece PAX Sling 150 cm light blue

2) Before use / after use

The stretcher is rolled out:



The belts are positioned outside the stretcher:

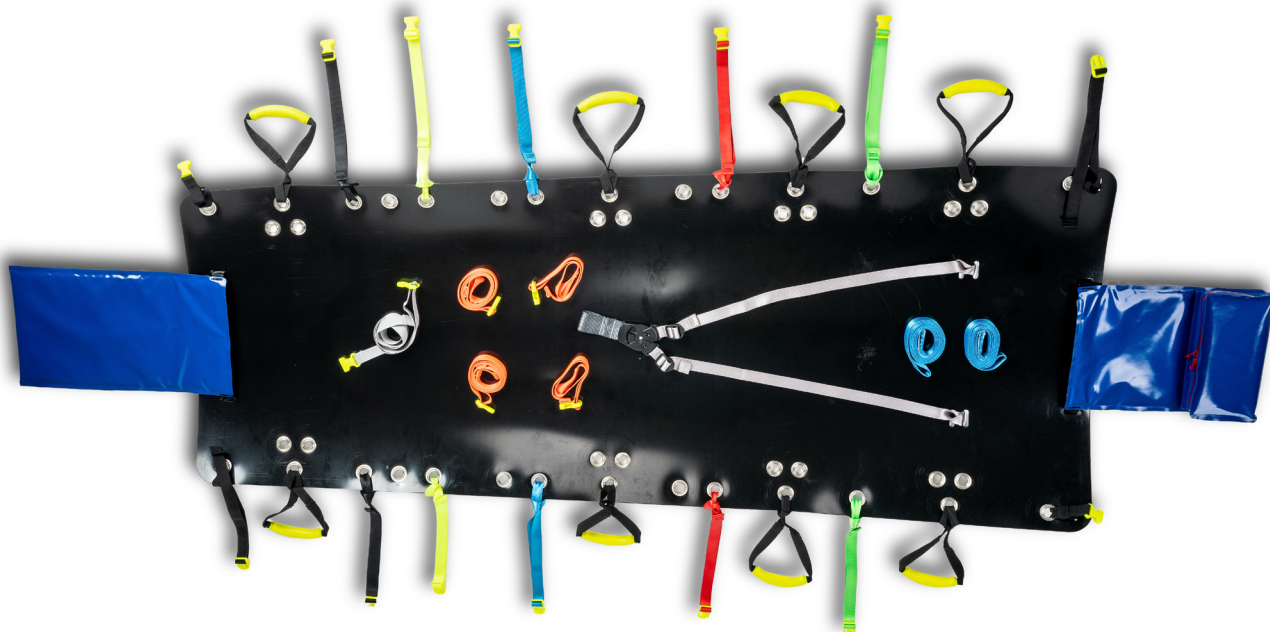



2) Before use / after use

The vacuum support is inserted into the stretcher:



PAX Vacuum Support for the LAS Stretcher



3) Packaging

The vacuum support bag and the tensioning strap are stored inside the transport bag of the stretcher:



The ends of the stretcher bag are brought together:



And secured by the integrated belt:





The additional slings and leg straps are stored in the zip pocket at the head end:





A) Preparing the vacuum support

Make sure that the PAX vacuum support is completely unfolded and not vacuumed.
Make sure that the vacuum valve is closed.

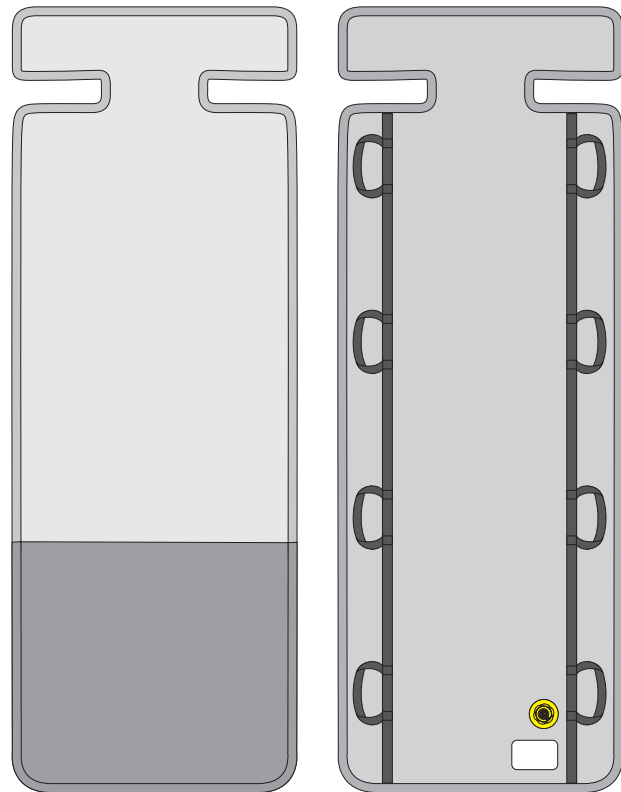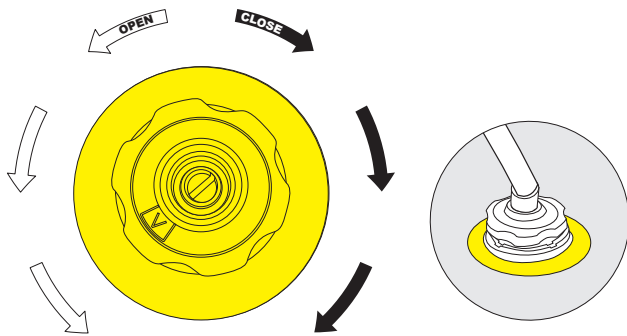

Front view

Rear view

B) Patient positioning

The PAX vacuum support enables improved patient positioning in the PAX Roll-Up Stretcher LAS.
The vacuum support is modelled to the patient to create a vacuum to promote a comfortable and safe position for the patient.
The vacuum support is formed around the patient to achieve a stable and comfortable position.

! Attention


The vacuum support must not be formed above/around the head or from below/around the feet, as this can lead to unwanted pressure on the spine.
Remove air and reshape the mattress.

The patient can now be positioned on the stretcher in a way that is suitable for transport. Close all buckles and tighten the fixing straps.
The straps are pulled taut in several steps from the foot end to the head end until the stretcher embraces the patient well and securely.

The transportation is carried out by at least 4 persons using the installed handles, 2 on each side.

If required, the stretcher can also be transported by means of the belt slings supplied.

For this purpose, the belt slings are looped through the left and right side of the head end of the stretcher and hung over the shoulder of the carrier.

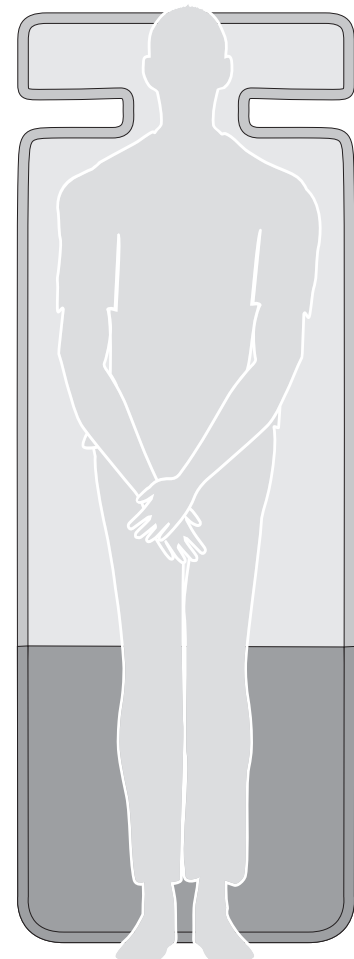

C) Preparing for transport

Mounting the optional leg loops (included in delivery) to adjust the vertical position of the patient in the stretcher.

The patient can be fixated in the correct vertical position on the stretcher by using the orange leg loops.

The short orange straps (6.0) are looped in through two eyelets in top handle reinforcement (red triangle with 3 eyelets) on the left and right side (position a).

The long orange straps (6.1) are looped through 2 eyelets in the handle reinforcement directly below (position b).

The ends of the straps with the buckle are then passed through the crotch above and below the patient and are then joined on the upper side of the patient.

The straps are pulled taut on both sides to adjust for correct patient position.




4) General information

Read these instructions before using the PAX Roll-Up Stretcher London Ambulance Service and ensure that they are fully understood.

The instructions must always be available in the respective national language.

The instructions must be made available to the user.

5) Information on the Label

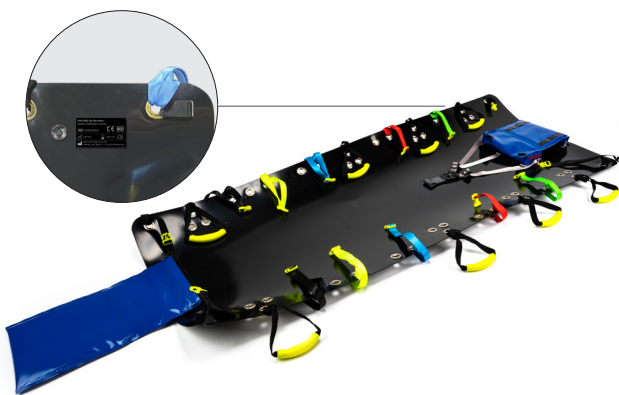

(A) Serial number

(B) Item number

(C) Article name

(D) Note icon (Read user notes!)

(E) Note icon (Article meets CE-standard)

(F) Note icon (Manufacturer
-name and -postal address)

6) Compatibility

The rescue system is made up of the individual components shown and may only be used with tested and approved components within the described operating conditions, as there **is a risk to life and limb if these are not complied with fully.**

7) Safety requirements

The instructions must be strictly observed for the protection of the user and the equipment!

The product marking must be completely legible!

It is essential that the user regularly checks all fastening and/or adjustment parts during use.

- Max carry load: 200 kg
- Application temperature -50°C to +50°C
- Store and transport in a dry place, protected from light
- Avoid contact with aggressive substances
- Observe danger from edges and rough surfaces
- Observe danger from electrical installations
- Observe danger due to moving machine parts
- Avoid heat and flames

8) Inspection and service life

A regular inspection according to country-specific requirements is necessary.

In Europe, an annual inspection by an expert is necessary.

Service life before the first use: 2 years with proper storage.

Service life from the first use: 5 years with proper storage.

The discard criteria depends on the product, its frequency of use and the external operating conditions.

Every piece of equipment loses its durability in the course of its service life. This is caused and accelerated by UV rays, heat, chemicals & aggressive atmospheres.

An inspection of the complete system and its individual components must be carried out at least once a year according to the regional requirements and carried out and documented by an expert.

A control card is available for this purpose under point 12.

9) Maintenance and storage

The PAX Roll-Up Stretcher London Ambulance Service should be checked regularly, preferably before each use, for damage and functionality.

In addition to the stretcher itself, also check the belt strap and handles.

If you discover any damage to these components, replace them immediately (if possible).

Check the vacuum support of the stretcher for leaks.

The vacuum support can be checked for leakage using the following method:

Cover the entire mattress with soapy water and blow a little air into it using a pump.

Do not inflate the mattress too hard! If bubbles appear in any part of the mattress, it is a leak that must be repaired. Use the PAX vacuum repair set for this.

If the damage is too large, please send the vacuum support to PAX for repair.

10) Cleaning

Clean the PAX Roll-Up Stretcher London Ambulance Service in lukewarm water (with neutral soap if necessary).

Rinse well afterwards.

Dry at room temperature.


Never place in a tumble dryer, in direct sunlight or near heaters!




Only
neutral soap




Do not wash
above 40°C

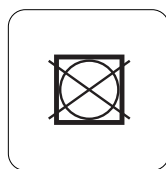

No
tumble dryer

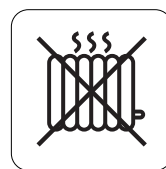

No
heaters

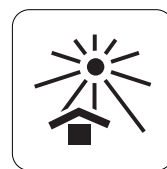

No
direct sunlight

11) Modification, repair

Changes and repairs may only be carried out by the manufacturer (does not apply to the attachment of spare parts).

If modifications or repairs are carried out independently on the PAX Roll-Up Stretcher London Ambulance Service, the operating licence will expire.



12) Control card

This control card must be filled in completely by the expert for PPE and medical devices during the annual expert inspection.

The control list does not claim to be complete and does not release the expert from assessing the overall condition.

Control card

Item name	
Serial number	
Date of manufacture	
Date of purchase	
Date of first use	
Product to be sampled on	

Year	Date	Name + Signature of the qualified person	Identified damage	Passed	Not passed	Measures
1						
2						
3						
4						
5						
6						
7						



X-CEN-TEK GmbH & Co. KG

Westerburger Weg 30
26203 Wardenburg
Deutschland / Germany

Tel. / Phone: +49 (0)44 07 - 7 14 76 - 0
Fax / Fax- : +49 (0)44 07 - 7 14 76 - 99
Email: info@pax-bags.de



www.pax-bags.com

