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PAX Infu-Warm-System

Insulating bag for storing intravenous drips in very cold environment, incl. heating gel pad which can be reused several times and insulated tube for the drip system.

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PAX Infu-Warm-System:

The isolation bag

The PAX Infu-Warm-System warms up an infusion during application on scene and works independent from electric sources. Keeping a patient warm is a big challenge during an operation. This challenge is increasing by sinking temperatures and during application of „cold" infusions what can lead to compications The PAX Infu-Warm-System offers the medication of a pre-warmed infusion. The warmth is generated by a reusable gel pad. The bag is lined with a warmth reflexive foil. The infusion system is insulated by a hose cover.

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PAX Infu-Warm-System

Filling level of the infusion can be observed by a see-through window without having to open the zipper. Only the velcro lid has to be opened for a short time. The bag itself remains closed. The PAX Infu-Warm-System is hung by a hand loop. It can be opened by a snap lock buckle.

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
The rear side of the bag is equipped with an additional compartment which contains an insulating hose for the infusion system included in the delivery volume.

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The delivery volume also includes a heating pad that can be reused several times. It is activated by only one movement of the hand. The heating pad is put into the mesh pocket determined to store it after it has been activated. Activation is done by a single kink of the metal plate and the exothermic reaction by which the liquid gel is turned into a crystalline material is started.

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The heating gel pad can be reused a hundredfold. (see instruction for use of the pad)

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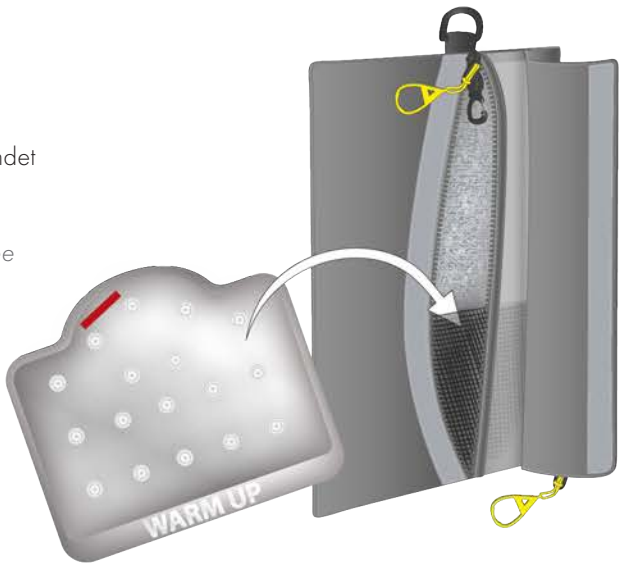The reuseable heating pad

It's ability for generating warmth at the exact temperature needed doing without external power sources where and when you need it is an important medium particularly with low temperatures during the operation. Only by pressing the plate the pads are heated up to approximately 50°C within seconds. If maintained correctly, they are reuseable a hundred times.

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Instructions for use:

Always put the heating pad on a flat ground before activation in order to ensure an uniform distribution of the liquid. Put the plate of the interior of the heating pad between thumbs and forefingers of both hands. Kink the plate once or several times (do not be afraid-it does not break!) and you will see how the liquid in the pad changes and you feel that it is heated up. Knead the pad after activation and complete crystallisation for a short time so that it can develop it's full heat and becomes comfortably soft. Now you can put the heating pad in the mesh pocket of the PAX Infu-Warm-System determined to store it. If the heating pad is used for something else, do not put the heating pad direct onto the skin.

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Regeneration of the heating pad for reuse:


1. *Boil some water in a pot (size to be chosen according to size of pad and quantity).*
2. *Put a thin towel in the pot an put the pad/s in the boiling water.*
 Boil the pad/s for approx. 10 - 15 minutes until the content is clear and liquid again.
3. *Take the pad out of the water. If there are still solid particles in the pad or if the liquid is getting solid again afer a short time, shake the pad a bit and put it into boiling water again for a few minutes.*
4. *If safe disinfection is required during use, the heating pad must be boiled for 30 minutes before reactivation.*
5. *Attention! Never use the heating pad directly after boiling since it is very hot then and with direct skin contact there is a danger of burn! After a short cool-down time your heating pad is ready to be used again.*

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Safety indication:

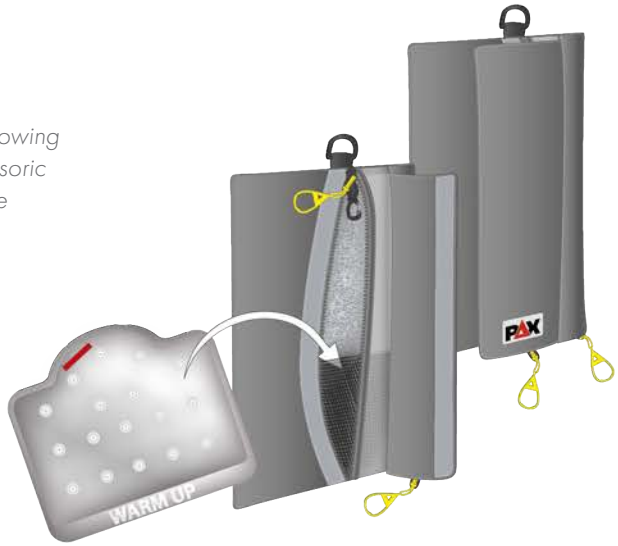If the heating pad is used beyond the PAX Infu-Warm-System, following advices are to be observed: Persons with blood circulation or sensoric problems should consult their medic prior to use. The content of the pad is not dangerous, not flammable, not toxic, not acidly and is ecologically compatible. Ingredients are approved according to TÜV (a German Technical Inspectation Authority)! If the pad is bursting by accident, the liquid becomes solid immediately. This can be easily removed with water. Soiled clothes or textiles can be normally washed. Unsuitable for consumption. No toy for children!

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The contents of this manual are protected by copyright. Any use, even in extracts, requires written permission.
X-CEN-TEK GmbH & Co. KG