




it's all about performance

PAX-Clip-In-System 3.0

Full concentration on administering medication.

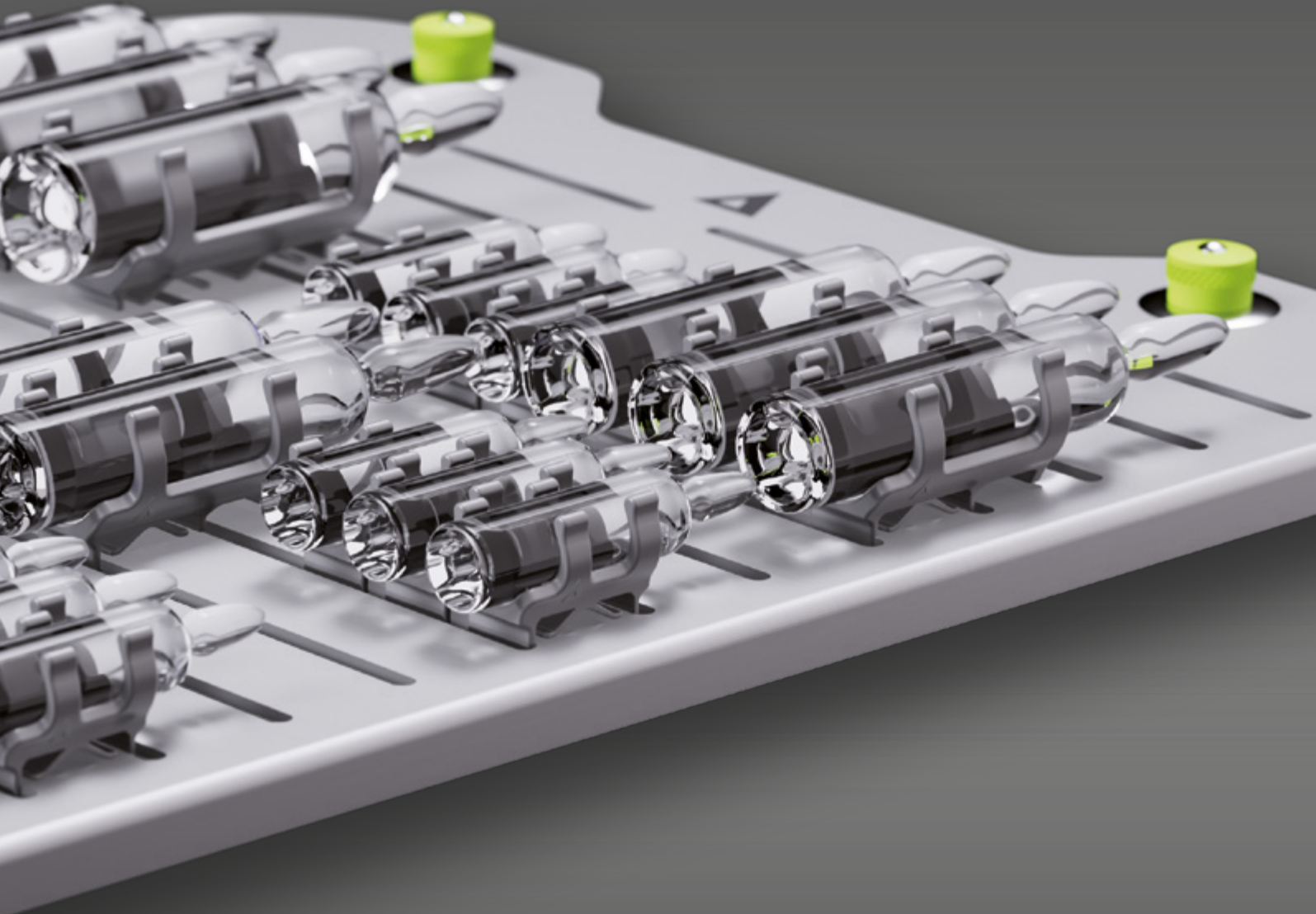

With the new PAX-Clip-In-System 3.0, you can noticeably improve processes. Efficiency and convenience combined in a sophisticated system.

▲ Increased patient and staff safety

▲ Simple installation

▲ Space-saving storage

▲ Increased clarity

▲ Reduced system weight

www.pax-bags.com





it's all about performance

Item no. 286274510

PAX ampoule holder XL PCI POS 3.0

The PAX ampoule holder XL PCI POS 3.0 is designed to fit the PAX Organising System (POS) and has the same external dimensions as the XL module. The PAX-Clip-In-System is mounted in the base and lid and is fitted with PAX ampoule clips as required. The centre walls, which can be folded out to the side, hold syringes, cannulas and other accessories and at the same time protect the ampoules from high external pressure. The PAX ampoule holder is fitted with fleece strips on the outside of the base to secure it in suitable backpacks. The number of ampoules that can be accommodated varies between approx. 70 and 200, depending on the size of the ampoules and the level of organisation required.

24x33x10 2,00

10
RIP
tec PRO

Item no. 286254510

PAX ampoule holder P5/11 M PCI 3.0

This PAX ampoule holder fits the P5/11 backpack system and is held in the backpack with magnets. It corresponds in size to the PAX See-Through module „MLH“ or two See-Through modules „SH“. The number of ampoules to be accommodated varies between approx. 60 and 160 ampoules, depending on the size of the ampoules and the desired number of ampoules.

24x30x10 2,05

10
RIP
tec PRO

Item no. 286264510

PAX ampoule holder P5/11 L PCI 3.0

In the modular system of the PAX P5/11 backpack, this PAX ampoule holder corresponds to the space of two PAX See-Through modules „L“ and is also held in position with magnets. This PAX ampoule holder currently offers the largest space for medication (up to approx. 200 ampoules). The centre walls can be folded out to the side to accommodate syringes, cannulas and other accessories, while also protecting the ampoules from external pressure.

30x36x10 2,60

10
RIP
tec PRO

Item no. 289194510

PAX function module medication PCI 3.0

Functional module for the safe storage of medication with a PAX-Clip-In-System carrier plate for FM medication integrated in the lid. The practical compartmentalisation in the body is equipped with transparent lids for clear storage of accessories. It is possible to attach an optionally available PAX handle for module P5/11 2.0 (Item no. 2788503..) on the outside.

28x36x10 1,26

10
RIP
tec PRO



it's all about performance

Item no. 289184510

PAX function module MLH medication PCI 3.0

Functional module for the safe storage of medication with a PAX-Clip-In-System carrier plate P5/11 M integrated in the lid. The practical compartmentalisation in the body, equipped with transparent lids, ensures clear storage of accessories. An optional PAX handle for module P5/11 2.0 (Item no. 2788503..) can be attached to the outside.

24x32x10 0,91

10

RIP PRO



PAX-Clip-In-System Module assembly

Item no.	Article name	Dimension*	**S 3.0	**M 3.0	**L 3.0	**XL 3.0	***Carrier Plate
286274510	PAX ampoule holder XL PCI POS 3.0	24x33x10	15	15	10	10	2x XL POS
286254510	PAX ampoule holder P5/11 M PCI 3.0	24x30x10	15	15	10	10	2x P5/11 M
286264510	PAX ampoule holder P5/11 L PCI 3.0	30x36x10	20	20	15	12	2x P5/11 L
289194510	PAX function module medication PCI 3.0	28x36x10	10	10	10	8	1x FM medication
289184510	PAX function module MLH medication PCI 3.0	24x32x10	10	10	5	6	1x P5/11 M

*in cm
**PAX-Clip-In-System ampoule clip
***PAX-Clip-In-System carrier plate

PAX-Clip-In-System Vehicle bracket

Item no.	Article name	Abmessung*	**S 3.0	**M 3.0	**L 3.0	**XL 3.0	***Trägerplatte
289384312	PAX-Clip-In-System Vehicle bracket P5/11 M	22,7x29	10	10	5	6	P5/11 M
289394312	PAX-Clip-In-System Vehicle bracket P5/11 L	28,5x35	10	10	10	8	P5/11 L
289404312	PAX-Clip-In-System Vehicle bracket XL POS 3.0	22,8x33	10	10	10	8	XL POS

*in cm
**PAX-Clip-In-System ampoule clip
***PAX-Clip-In-System carrier plate


The clip makes the difference

The new design of the PAX-Clip-In-System ampoule clips keeps everything organised and saves space. Four differently-sized clips ensure a secure hold, easy assembly and a clearly organised arrangement.


PAX-Clip-In-System ampoule clip S 3.0




PAX-Clip-In-System ampoule clip M 3.0

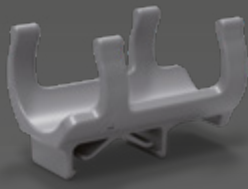

PAX-Clip-In-System ampoule clip L 3.0

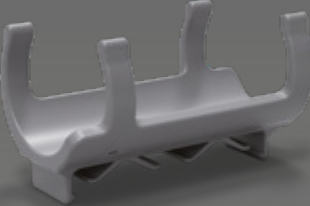

PAX-Clip-In-System ampoule clip XL 3.0

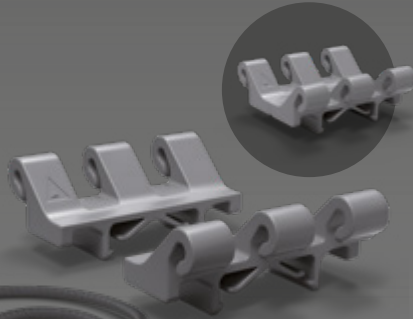

PAX-Clip-In-System carrier plate P5/11 M



PAX-Clip-In-System carrier plate P5/11 L



PAX-Clip-In-System carrier plate XL



PAX-Clip-In-System carrier plate FM medication



PAX-Clip-In-System ampoule clips 3.0


		*S 3.0	*M 3.0	*L 3.0	*XL 3.0
Glass break ampoules	1 ml	x			
	2 ml	x			
	5 ml		x		
	10 ml			x	
Ampoules					x
Plastic ampoules**			x	x	x
Plastic vials		x			

*PAX-Clip-In-System ampoule clip
**Depending on manufacturer

PAX-Clip-In-System ampoule clip sets

Item no.	Article name	Number of clips
289340010	PAX-Clip-In-System ampoule clip S 3.0 set	Set of 20 clips S 3.0
289350010	PAX-Clip-In-System ampoule clip M 3.0 set	Set of 20 clips M 3.0
289360010	PAX-Clip-In-System ampoule clip L 3.0 set	Set of 20 clips L 3.0
289370010	PAX-Clip-In-System ampoule clip XL 3.0 set	Set of 20 clips XL 3.0*

*Plus 1 m elastic cord 1 mm
(10 total clips)

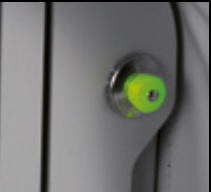




Further
information:



X-CEN-TEK GmbH & Co. KG
Westerburger Weg 30
26203 Wardenburg / Germany

+49 (0)44 07 - 7 14 76 - 0
info@pax-bags.de

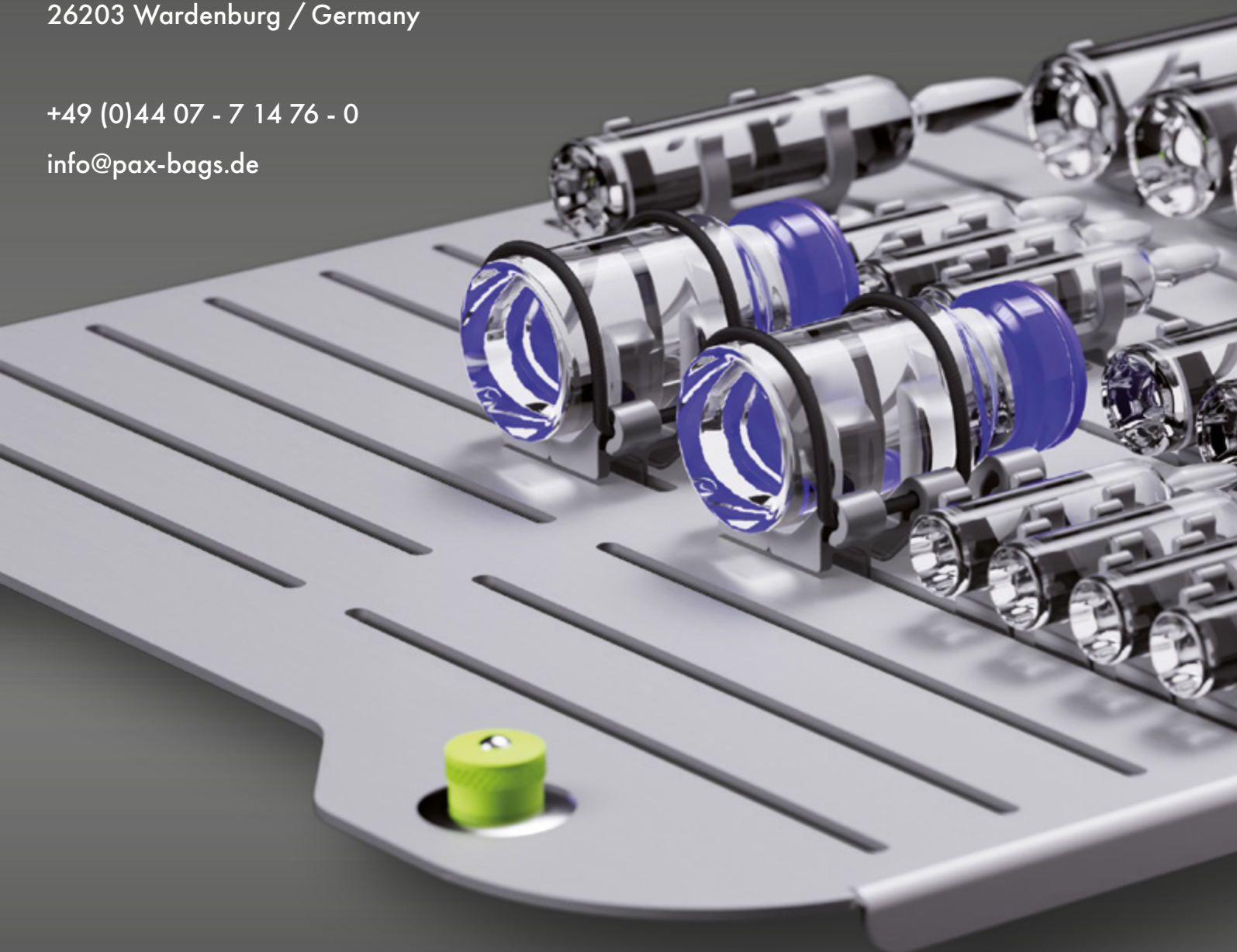

Errors or changes for technical reasons
that serve to improve the product cannot be ruled out.
All products are delivered without content.